

Water Quality Monitoring: Total Dissolved Oxygen Assessment in Different Water Samples by Winkler's Method

1. Importance of Dissolved Oxygen (DO): Knowledge of DO concentration in seawater is often necessary in environmental and marine sciences. It is used by oceanographers to study water masses in the ocean. It provides the marine biologist a means to measure primary production, particularly in laboratory cultures. For the marine chemist, it provides a measure of the redox potential of the water column. DO is also an important factor in corrosion. Oxygen is poorly soluble in water. The solubility of oxygen decreases with increase in concentration of the salt and hence, solubility of DO is lesser in saline water. The amount of DO at 100% saturation at sea level is 9.03 mg/L (at 20 °C) and is sufficient to sustain aquatic life. DO is used as an indicator of the water body's health, where higher DO concentrations are correlated with high productivity and little pollution. Dissolved oxygen is usually determined by Winkler's method.

2. How does the Winkler Method work? Winkler's Method uses titration to determine DO in water sample. A sample bottle is filled completely with water (no air is left to skew the results). DO in the sample is then "fixed" by adding a series of reagents that form an acid compound that is then titrated with a neutralizing compound that results in a colour change. The point of colour change is called the "end-point," which coincides with the DO concentration in the sample. DO analysis is best done in the field, as the sample will be less altered by atmospheric equilibration.

3. Applications:


DO analysis can be used to determine the health or cleanliness of a lake or stream, amount and type of biomass a freshwater, the amount of DO that a system can support, and the amount of decomposition occurring in the lake or stream.

Expt. No.: 2

Date: 2021/02/04

Experiment	Water Quality Monitoring: Total Dissolved Oxygen Assessment in Different Water Samples by Winkler's Method
Problem definition	Dissolved oxygen (DO) is essential to living organisms in water but harmful if present in boiler feed water leading to boiler corrosion.
Methodology	Winkler's titration method is used to assess DO in water.
Solution	Estimation of total dissolved oxygen in different water samples.
Student learning outcomes	Students will learn to a) perform Winkler's titration method b) assess the total dissolved oxygen in different water samples

Principle: Estimation of dissolved oxygen (DO) in water is useful in studying corrosion effect of boiler feed water and in studying water pollution. DO is usually determined by Winkler's titration method, which is based on the fact that DO oxidize potassium iodide (KI) to iodine (I₂). Since DO in water is in molecular state, as such it cannot oxidize KI. Hence, manganese hydroxide is used as an oxygen carrier to bring about the reaction between KI and Oxygen. Manganese hydroxide, in turn, is obtained by the action of NaOH on MnSO₄.

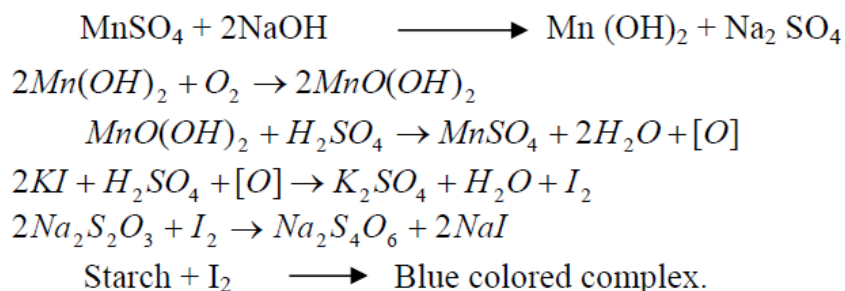

The liberated iodine (I₂) is titrated against standard sodium thiosulphate (Na₂S₂O₃) solution using starch as indicator.

Requirements:

Reagents and solutions: Standard buffer of pH 7, standard potassium dichromate (0.01 N), sodium thiosulphate solution, 10% KI solution, alkali KI solution (KI + KOH in water), conc. H₂SO₄, manganese sulphate, starch solution as indicator.

Apparatus: Conical flask, Burette, Measuring flask, Beakers.

Procedure:

Titration 1: Standardization of Sodium Thiosulphate

Rinse and fill the burette with given sodium thiosulphate solution (Bottle B). Pipette out 20 mL of 0.01N $K_2Cr_2O_7$ solution (Bottle A) into a clean conical flask. To this, add 5 mL H_2SO_4 (1/2 T.T.), 5 ml (1/2 T.T.) of 10% KI, and titrate against sodium thiosulphate solution. When the solution becomes straw yellow colour, add starch indicator and continue the titration. End point is the disappearance of bluish brown colour. Repeat the titration to get concordant value.

Titration 2: Estimation of Dissolved Oxygen

Using a measuring cylinder, add 100 mL of sample water in a conical flask. Further, add 2 mL of $MnSO_4$ and 2 mL of alkali KI solution and shake well for the rough mixing of the reagents. Set aside the flask for few minutes to allow the precipitate to settle down, and then add 2 mL of conc. H_2SO_4 for complete dissolution of the precipitate. Then, titrate against std. sodium thiosulphate solution. When the solution turn into light yellow, add starch indicator. End point is the disappearance of bluish brown colour. Repeat the titration to get the concordant value. Calculate the strength of dissolved oxygen from the titer value. Based on that, calculate the amount of DO in the given water sample.

OBSERVATION AND CALCULATIONS

Titration - I: Standardization of Sodium Thiosulphate

S. No.	Volume of $K_2Cr_2O_7$ (V_1 , mL)	Burette reading (mL)		Volume of sodium thiosulphate (V_2 , mL)
		Initial	Final	
1.	20 ml	0	16.6	16.6
2.	20 ml	0	16.7	16.7
3.		0	16.7	16.7
Concordant value				16.7

Calculations:

Volume of potassium dichromate $V_1 = 20$ mL

Strength of potassium dichromate $N_1 = 0.01$ N

Volume of sodium thiosulphate $V_2 = \dots 16.7 \dots$ mL (From Titration – 1)

Strength of sodium thiosulphate $N_2 = \dots \dots \dots ?$

$$V_1 N_1 = V_2 N_2$$

$$\therefore N_2 = V_1 N_1 / V_2$$

$$\begin{aligned} \text{Strength of sodium thiosulphate} = N_2 &= 20 \times 0.01 / V_2 = \dots \dots \dots 0.01197 \text{ N} \\ &\approx 0.012 \text{ N} \end{aligned}$$

Titration – II: Estimation of Dissolved Oxygen

S. No.	Volume of water sample (V ₁ , mL)	Burette reading (mL)		Volume of sodium thiosulphate (V ₂ , mL)
		Initial	Final	
1.	100 ml	0	7.2	7.2
2.	100 ml	0	7.2	7.2
3.				
Concordant value				7.2

Calculation:

Volume of sodium thiosulphate V₂ = 7.2 mL (From Titration – 2)

Strength of sodium thiosulphate N₂ = 0.00976 N (From Titration – 1 calculation)

Volume of water sample taken V₁ = 100 mL

Strength of given water sample N₁ = ?

$$V_1 N_1 = V_2 N_2$$

$$N_1 = V_2 \times N_2 / 100$$

$$= 0.000862 \text{ N}$$

Amount of DO (ppm) = normality × equivalent weight of O₂ × 1000 mg/L of the given water sample.

$$= 8.62 \times 10^{-4} \text{ N} \times 8 \times 1000 \text{ mg/L}$$

$$= 6.89 \text{ ppm.}$$

Result: Amount of DO in the given water sample = 6.89 ppm.